
Structural MRI Unwarping Using CMTK¹

Release 1.4

Torsten Rohlfing

May 30, 2013

Neuroscience Program, SRI International, Menlo Park, CA

Abstract

This document describes the workflow for unwarping structural MR images, in particular T_1 -weighted SPGR and MP-RAGE scans, using reference scans of the Magphan[®] EMR051 Quantitative Imaging Phantom (a.k.a. ADNI Phantom) and the tools of the Computational Morphometry Toolkit (CMTK).

Contents

1	Introduction	2
2	What You Need and Where You Get It	2
2.1	Structural Imaging Phantom	2
2.2	CMTK	3
3	Step-by-Step	3
3.1	Imaging	3
3.2	DICOM Image Stacking	3
3.3	Phantom Detection	5
3.4	Fault-Tolerant Phantom Detection	5
3.5	Create Unwarping Transformation	7
3.6	Correction With and Without Unwarping	7
4	Concluding Remarks	8
A	Example Data	10
B	Example XML Phantom File	10
C	Phantom Geometry	13

¹This document is licensed under the Creative Commons Attribution License Version 3.0.

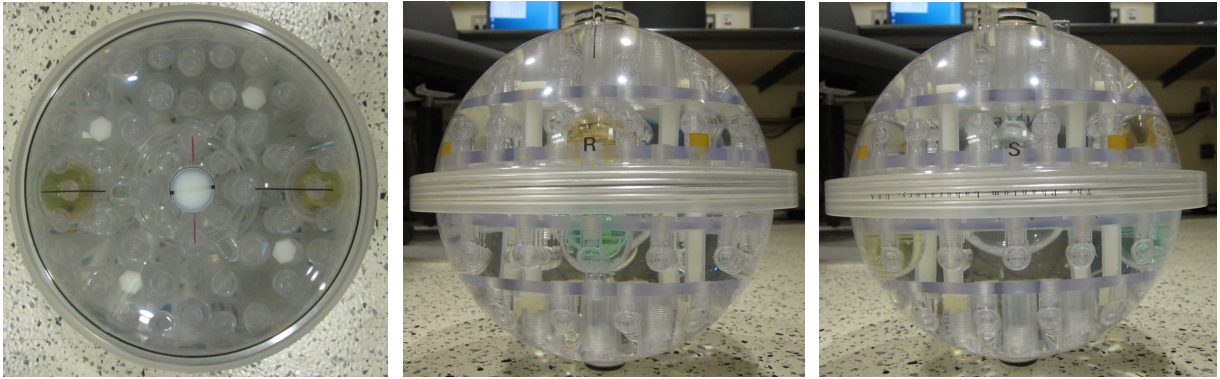


Figure 1: The Magphan® EMR051 Quantitative Imaging Phantom, seen from “anterior,” “right,” and “superior” direction (in “patient coordinates,” pictures from left to right).

1 Introduction

Nonlinear distortions and drift of scale calibration can severely confound imaging-based studies, especially when run on multiple scanners or even at multiple imaging sites.

The Magphan® EMR051 Quantitative Imaging Phantom (also known as the “ADNI Phantom” after the Alzheimer’s Disease Neuroimaging Initiative) provides a technical solution for detecting, and ultimately correcting, the effects of nonlinearity as well as drift over time, whether on a single scanner or across devices.

Unfortunately, software tools to make practical use of the phantom are not freely available [1] (see also [3]).

We describe in this article for the first time a workflow implementing the correction of scanner miscalibration and nonlinearities *using only freely available data and software tools*. Example image data is provided with the article. Source code for all software tools is available from <http://nitrc.org/projects/cmtk/>.

2 What You Need and Where You Get It

We will assume you already have an MR scanner. If not, pick up a used one on eBay. Make sure you keep an eye on the shipping charges.

2.1 Structural Imaging Phantom

You will need the “Magphan® EMR051 Quantitative Imaging Phantom” (also known as “The ADNI Phantom” [1]), http://www.phantomlab.com/products/magphan_adni.php. Photos of one such phantom are shown in Fig. 1.

The phantom can be purchased from the manufacturer, The Phantom Laboratory, P.O. Box 511, Salem, NY, 12865-0511 USA. These are expensive – hope you saved some cash when you bought your scanner.

2.2 CMTK

Unlike the previous two items, The Computational Morphometry Toolkit (CMTK) is **free**, and that's as in both free beer and free speech. CMTK is available both in source code, licensed under the GPLv3, and as pre-compiled binary distributions from <http://nitrc.org/projects/cmtk/>. If you are using NeuroDebian, you can also install CMTK directly.

We shall assume that CMTK has been installed such that its tools can be run as

```
cmtk <tool> <arg1> <arg2> ...
```

You will need CMTK release 2.2.0 or later. Earlier versions do not support phantom detection of landmark-based nonlinear deformations.

3 Step-by-Step

3.1 Imaging

The imaging phantom should be placed in the scanner according to manufacturer instructions and T_1 -weighted images should be acquired using the site-preferred imaging sequence (e.g., SPGR or MP-RAGE).

It is important to make sure that the entire phantom is contained within the acquisition field of view. Add slices if necessary.

Example DICOM images of the phantom acquired with “sagittal” slice orientation and full coverage are shown in Fig. 2.

3.2 DICOM Image Stacking

Assume that the DICOM files containing the phantom images are stored in the “dicom/” directory. These are stacked into a 3D image in NIFTI-format using the following CMTK command:

```
cmtk dcm2image -O phantom.nii dicom/
```

Make sure only the DICOM files of the actual structural scans are in the “dicom/” directory, that is, not additional files such as scout images¹.

While CMTK supports a number of file formats, **only NIFTI or NRRD format should be used for storing the phantom image.** The same is true for any patient images to be unwarped. The reason for this limitation is that, of the supported formats, only NIFTI and NRRD preserve the physical scanner coordinates at which the images were acquired. This information is absolutely necessary to determine the correct spatial relationship between phantom and subject images. Without establishing this relationship, images cannot be unwarped.

¹Otherwise, you can change the output file name to “phantom%n.nii”, which will result in multiple, numbered output images. You will then need to identify the correct one.


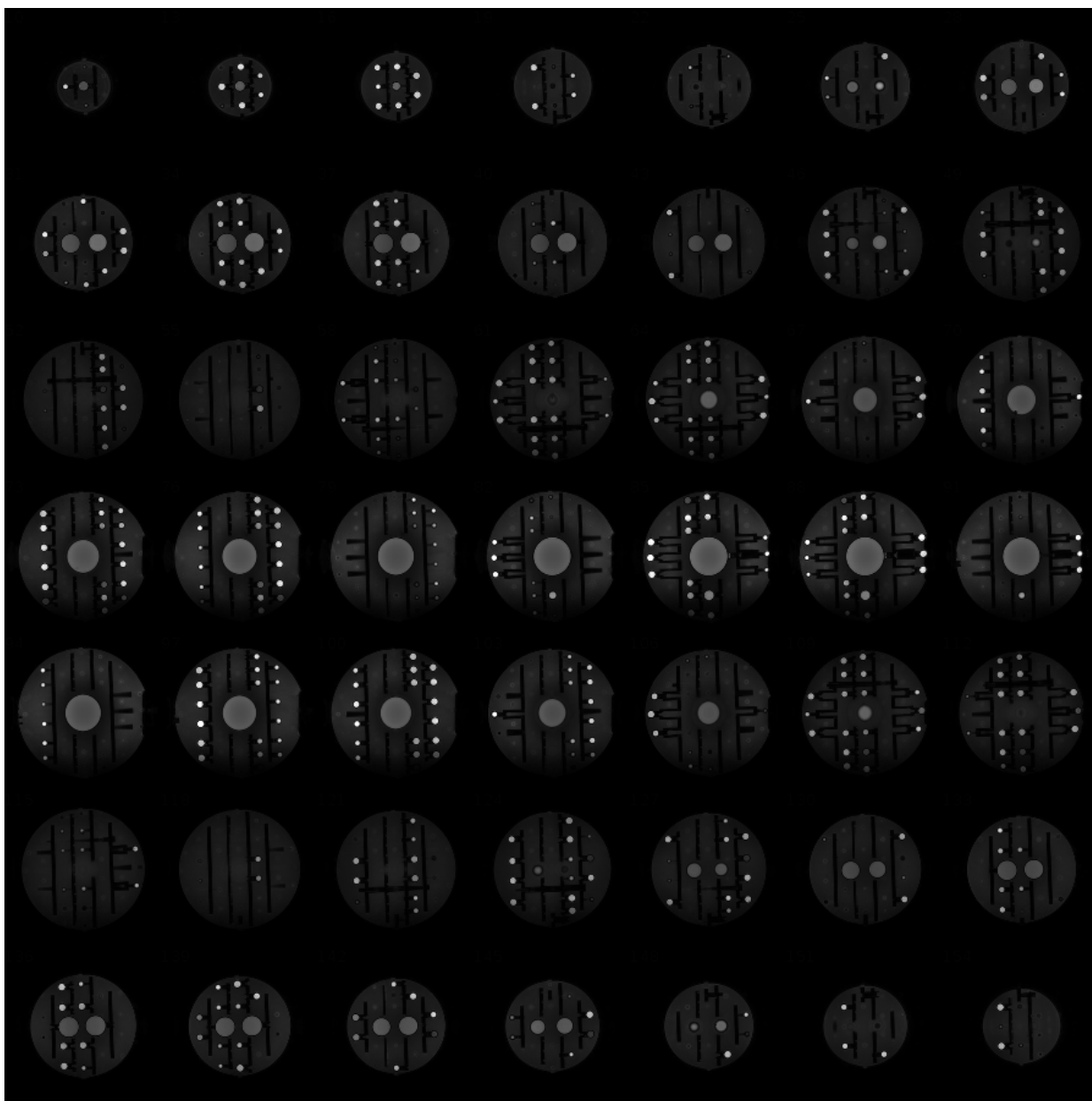


Figure 2: Thumbnails of selected DICOM images acquired using SPGR acquisition on a GE MRI scanner (total number of slices in the acquisition was 176).

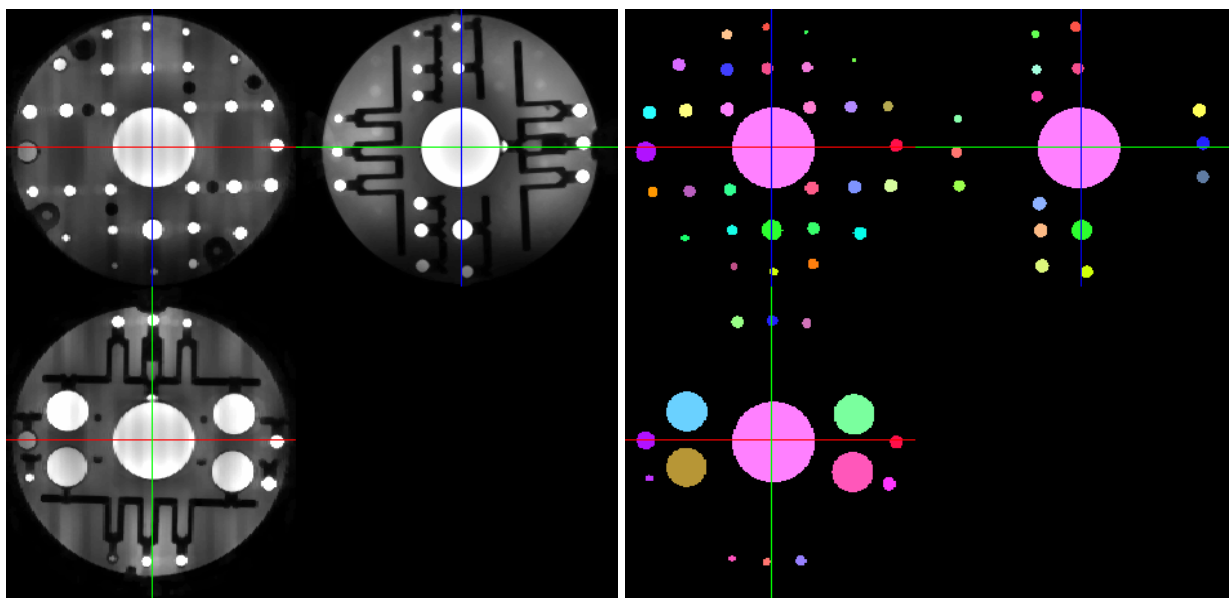

3.3 Phantom Detection

```
cmtk detect adni phantom phantom.nii phantom.xml
```

To verify landmark locations, the phantom detection tool can optionally write a label image matching the input image, in which the volume of each detected phantom sphere is marked with a unique label. To this end, add “`--write-labels labels.nii`” to the above command line to produce the labels file. See Fig. 3 for an example of three orthogonal slices from a phantom image with corresponding labeled spheres.

Not all MR scans of the imaging phantom may be of ideal quality, and sub-standard scans may cause phantom detection to fail. Examples are images with insufficient field-of-view, where not all phantom spheres are actually fully covered. Also the phantom itself may be degraded, resulting in insufficient contrast of some spheres.

By default, the phantom detection tool will terminate operation with an error message if problems such as missing or undetectable spheres are encountered. To produce potentially usable results in these cases,

however, the tool offers an optional “tolerant” mode of operation (beginning with CMTK Release 2.3.0). This mode is invoked by adding the “--tolerant” option to the command line:

```
cmtk detect_adni_phantom --tolerant phantom.nii phantom.xml
```

Note that in this mode, the tool can even potentially recover from otherwise fatal conditions, such as a missing 15 mm-diameter sphere (which is required by the standard phantom detection algorithm to establish gross phantom orientation). Instead, in tolerant mode, the detection tool will attempt to bootstrap orientation from the positions of the 30 mm-diameter CNR spheres, and the only requirement here is that the Orange CNR sphere appears with the highest intensity of all four CNR spheres (thus making detection robust even when doping of the CNR spheres is compromised with respect to one another).

In more technical detail, the phantom detection tool identifies the phantom orientation using the following procedure:

1. The 60 mm SNR sphere and all four 30 mm CNR spheres are detected independently using an FFT-based matched, bipolar filter.
2. The two 15 mm spheres are detected, again using a matched filter. Search for each sphere is constrained to a ring of appropriate radius around the center of the 60 mm sphere.
3. The angle is computed between the two lines from the SNR sphere to each of the 15 mm spheres. The same is done for the lines from the centroid of the four CNR sphere centers to each of the 15 mm spheres.
4. If the resulting angle is closer to 90 degrees for the CNR-spheres centroid, then the latter is used as the phantom center rather than the center of the SNR sphere (this makes the procedure robust against phantoms with broken-off SNR spheres).
5. If the angle between lines from the the selected phantom center to the two 15 mm spheres is within the range of 85 to 95 degrees, then the coordinate system defined by the 15 mm spheres and the phantom center is used to establish gross phantom orientation (via a least-squares rigid alignment).
6. If said angle is outside the valid range, this indicates that one or both of the 15 mm spheres are missing or could not be detected for other reasons. In this case, the brightest CNR sphere (“Orange”) is identified by comparing mean intensities across all four CNR spheres. The CNR sphere closest to is identified as the “Green” sphere. Together with the center of the SNR sphere, these two landmarks define a coordinate system that is used as the third and final alternative for establishing phantom orientation.

In summary, phantom detection as implemented above can tolerate all three of the following phantom defects: a) one or two missing 15 mm spheres, b) misplaced SNR sphere, c) out-of-order CNR sphere intensities (so long as the “Orange” sphere remains the brightest).

In order to support this level of robustness, the phantom detection tool does have to make one assumption about the way that the phantom is scanned, namely that the phantom was scanned with the correct side facing up (the A direction in patient coordinates). If the phantom was scanned upside-down, then detection will fail, unless the “--any-orientation” command line option is invoked (which, in turn, will lead to detection failures with defective phantoms).

3.5 Create Unwarping Transformation

From the phantom description file, we can create an unwarping transformation using the following CMTK command (to be entered on a single command line):

```
cmtk unwarp_image_phantom --final-cp-spacing 40 --levels 2 \  
phantom.xml phantom.nii unwarp.xform
```

This will generate a nonrigid transformation, “phantom_warp.xform,” in the space of the image “phantom.nii.” If a subject image, “subject.nii,” is to be unwrapped, the command should be changed to (again entered on a single command line):

```
cmtk unwarp_image_phantom --final-cp-spacing 40 --levels 2 \  
phantom.xml subject.nii unwarp.xform
```

Note that a separate transformation must be created for every image to be unwrapped, as each image has different physical coordinates.

For now, the only supported representation for the unwarping transformation is a B-spline free-form deformation [4], which requires specification of a final control point spacing (here: 40 mm) and is improved by specifying also a number of multi-resolution levels for the spline fitting process (here: 2 levels) [2]. Future versions of CMTK will implement a thin-plate spline, in which each landmark will be used directly as a control point.

3.6 Correction With and Without Unwarping

Using the deformation field fitted to the phantom landmark locations, the phantom image can be unwrapped as follows (entering the command on a single command line):

```
cmtk reformatx --sinc-cosine -o unwrapped.nii --floating phantom.nii \  
phantom.nii unwarp.xform
```

The phantom image, “phantom.nii,” is listed twice because it is both the floating (moving) and the reference (fixed) image. Again, if a separate subject image is to be unwrapped, the command changes to

```
cmtk reformatx --sinc-cosine -o unwrapped.nii --floating subject.nii \  
subject.nii unwarp.xform
```

But do you really want to unwrap the image in the first place?

Here’s something to consider — unwarping the image, phantom or subject, requires interpolation and thus introduces smoothing, even though we are using a sinc-kernel above. That’s fine if you really need to know what anatomy is exactly where, say to use the unwrapped image for guiding implantation of a deep-brain stimulator. **(Except, it is not fine, because CMTK is not FDA approved, so don’t do it!)**

But what if you only want to unwrap the image to get the correct tissue or region volumes somewhere down your processing pipeline? Then you might actually be making things worse, for example for tissue segmentation, by blurring the image ever so slightly.

Instead, what you really want to know is how the volume of each pixel changes due to distortion so you can compute correct region and tissue volumes. Well, we can do that without interpolation.

What we need for this correct is a map that tells us, for each pixel, how much larger (or smaller) its volume really is, relative to the ideal volume as prescribed in the imaging protocol. Mathematically speaking, we need the map of Jacobian determinants of the inverse unwarping transformation in the coordinate space of the acquired image. Unfortunately, we also need this map deformed using the inverse of the transformation itself.

CMTK can give us this exact type of map using the following command (here, for a subject image, but analogous for the phantom image itself if you care):

```
cmtk reformatx -o pxvolume.nii subject.nii --inverse unwarp.xform --jacobian \
--inverse unwarp.xform
```

Equivalently, we can use the following two commands to avoid at least the second numerical inversion of the nonrigid transformation:

```
cmtk reformatx -o temporary.nii subject.nii --inverse unwarp.xform --jacobian \
unwarp.xform
cmtk imagemath --in temporary.nii --one-over --out pxvolume.nii
```

This makes use of the fact that the Jacobian determinant of forward and inverse transformations are related by

$$J_{T^{-1}}(T(\vec{x})) = 1/J_T(\vec{x})$$

for a transformation T and any location \vec{x} . The former command above computes the left-hand side of this equation, whereas the latter two-command sequence computes its right-hand side.

The best way, however, to create a Jacobian volume correction map is to avoid inversion of the unwarping transformation altogether by simply fitting a deformation to the “inverse” phantom landmarks instead, i.e., to the landmarks with detected and expected positions exchanged. The fitting tool provides the “--fit-inverse” option for just this purpose:

```
cmtk unwarp_image_phantom --fit-inverse --final-cp-spacing 40 --levels 2 \
phantom.xml phantom.nii unwarp_inverse.xform
```

Then we can simply evaluate the Jacobian of the fitted inverse transformation to obtain the volume correction map:

```
cmtk reformatx -o pxvolume.nii subject.nii --jacobian unwarp_inverse.xform
```

Either way, the resulting image, “pxvolume.nii,” contains at each pixel the true volume of that pixel. By multiplying these values with tissue volumes resulting from segmentation of the original image, we can obtain distortion-corrected volumes *without actually unwarping and thus blurring the original image*.

4 Concluding Remarks

At the time of writing, several of the software tools described herein (CMTK’s “detect_adni_phantom” and “unwarp_phantom_image” tools) are still considered work-in-progress. We are making these tools

and this article available to provide the community with an opportunity to evaluate our software in a timely fashion. Production use, however, it not recommended at this time.

Acknowledgments

This work was supported by the National Institute on Alcohol Abuse and Alcoholism (NIAAA) under Grant No. U01 AA021697, National Consortium for Alcohol and NeuroDevelopment in Adolescence (N-CANDA), Data Analysis Component.

Despite the name “ADNI Phantom” none of the materials used in this article and none of the phantom-related code in CMTK make use of any ADNI data or documentation, other than a cursory inspection of Ref. [1]. In particular, the geometric specifications of the phantom itself were derived from the phantom manual, available from http://www.phantomlab.com/library/pdf/magphan_adni_manual.pdf The table representing the geometry in CMTK’s source code is included in Appendix C and can be re-used under the terms of the CC-BY-3.0 license.

Use of the materials and software tools described in this article does not establish a requirement to acknowledge ADNI on the author list of subsequent publications [3].

References

- [1] J. L. Gunter, M. A. Bernstein, B. J. Borowski, C. P. Ward, P. J. Britson, J. P. Felmlee, N. Schuff, M. Weiner, and C. R. Jack. “Measurement of MRI scanner performance with the ADNI phantom.” *Medical Physics*, **36**(6):2193–2205, 2009. <http://dx.doi.org/10.1118/1.3116776>. 1, 2.1, 4
- [2] S. Lee, G. Wolberg, and S. Y. Shin. “Scattered data interpolation with multilevel B-splines.” *IEEE Transactions on Visualization and Computer Graphics*, **3**(3):228–244, 1997. <http://dx.doi.org/10.1109/2945.620490>. 3.5
- [3] T. Rohlfing and J.-B. Poline. “Why shared data should not be acknowledged on the author byline.” *NeuroImage*, **59**(4):4189–4195, 2012. <http://dx.doi.org/10.1016/j.neuroimage.2011.09.080>. PMID 22008368. 1, 4
- [4] D. Rueckert, L. I. Sonoda, C. Hayes, D. L. G. Hill, M. O. Leach, and D. J. Hawkes. “Nonrigid registration using free-form deformations: Application to breast MR images.” *IEEE Transactions on Medical Imaging*, **18**(8):712–721, 1999. <http://dx.doi.org/10.1109/42.796284>. 3.5

A Example Data

This article is provided with an example image (“phantom.nii.gz”) of an actual phantom. Running CMTK’s phantom detection tool on this image (with increased verbosity level for some interesting statistics) yields:

```
> cmtk detect_adni_phantom --verbose-level 2 --write-labels labels.nii \
    phantom.nii.gz phantom.xml
INFO: landmark fitting error average = 0.725241 maximum = 1.77297
      maxErrName = 10mm_0_12 maxErrLabel = 19
INFO: detected and matched 160 out of 160 expected landmarks.
```

This produces the phantom description file, “phantom.xml.” The contents of this file are explained below in [Appendix B](#).

We also see that all landmarks were successfully detected and the average linear transformation fitting residual over all landmarks was 0.72 mm. The sphere with the maximum residual is “10mm_0_12,” i.e., the 12th of the 10 mm spheres in Plane 0 of the phantom. (This sphere is labeled as ROI #19 in the phantom label file, “labels.nii.”

Using the XML phantom description, we can generate a nonrigid transformation for phantom unwarping:

```
> cmtk unwarped_image_phantom --final-cp-spacing 40 --levels 2 phantom.xml \
    phantom.nii.gz phantom.xform
```

Finally, we can use this transformation, “phantom.xform,” to create the unwarped phantom image:

```
cmtk reformatx -o unwarped.nii --floating phantom.nii.gz phantom.nii.gz \
    phantom.xform
```

B Example XML Phantom File

The `detect_adni_phantom` tool creates an XML file that contains a description of the phantom detected in a given image. An example of the contents of this file is shown below.

First, the file contains the XML header and the name of the represented phantom:

```
<?xml version="1.0" encoding="utf-8"?>
<phantom>
<phantomType>MagphanEMR051</phantomType>
```

Next, the file contains the estimated signal-to-noise ratio (based on the phantom’s SNR sphere) and four estimates of contrast-to-noise ratio (each based on one of the four CNR spheres and its contrast relative to the SNR sphere):

```
<snr>18.304893</snr>
<cnr>14.469221 35.670155 29.324328 28.170110</cnr>
```

Following is the list of detected landmarks – all coordinates are given in physical image coordinates, derived ultimately from the DICOM headers but represented in “RAS” coordinates (which means that x and y coordinates are negative relative to DICOM’s “LPS” coordinates):

```
<landmarkList coordinates="physical" space="RAS" count="165">
```

The `count` attribute provides the number of successfully detected landmarks (which may be lower than 165 in “tolerant” phantom detection mode).

Immediately following, for each landmark, its unique name, expected location, detected location, precision flag, and fitting residual are stored. The landmark name is essentially arbitrary and based on the table of phantom landmarks as shown in Appendix C.

The *expected location* is where the sphere center should be located in the phantom image if the transformation between ideal phantom and image was rigid (i.e., no scale, no shear, no nonlinearity).

The *detected location* is where the sphere center was actually detected in the image.

The *precision flag* specifies whether the ideal location of this sphere should be considered as precise based on the phantom construction (e.g., the SNR and CNR spheres have manufacturing tolerances that make their locations unprecise). Only landmarks with this flag set to “yes” should be used for registration.

The *fitting residual* is the Euclidean distance between the detected landmark location and the expected location according to a *linear fit* of all landmarks. Note that here, the linear fit may include anisotropic scale and shears, i.e., the residual is *not* simply the distance between expected and detected location as stored in the file. The purpose of this residual is to allow the transformation fitting tool to exclude outliers based on a threshold of allowable residual.

```
<landmark>
  <name>SNR</name>
  <expected>-2.081532 29.672649 -4.689435</expected>
  <detected>-4.846132 28.367746 -4.338875</detected>
  <isPrecise>no</isPrecise>
  <fitResidual>3.078036</fitResidual>
</landmark>
<landmark>
  <name>15mm@90mm</name>
  <expected>86.924671 28.899446 -7.279399</expected>
  <detected>87.598368 29.497402 -7.207202</detected>
  <isPrecise>no</isPrecise>
  <fitResidual>1.000909</fitResidual>
</landmark>
<landmark>
  <name>15mm@60mm</name>
  <expected>-3.739912 29.019053 -64.733024</expected>
  <detected>-3.855809 30.500179 -65.321349</detected>
  <isPrecise>no</isPrecise>
  <fitResidual>1.687546</fitResidual>
</landmark>
<landmark>
  <name>CNR-Orange</name>
```

```

        <expected>57.536139 46.089755 -5.790284</expected>
        <detected>58.235312 51.039218 -6.136547</detected>
        <isPrecise>no</isPrecise>
        <fitResidual>4.948012</fitResidual>
    </landmark>
    <landmark>
        <name>CNR-Red</name>
        <expected>-62.383089 43.226012 -2.446989</expected>
        <detected>-63.540399 48.950964 -2.568987</detected>
        <isPrecise>no</isPrecise>
        <fitResidual>5.653304</fitResidual>
    </landmark>
    <landmark>
        <name>CNR-Yellow</name>
        <expected>-61.699203 13.255542 -3.588586</expected>
        <detected>-62.273718 6.096574 -3.753603</detected>
        <isPrecise>no</isPrecise>
        <fitResidual>7.112524</fitResidual>
    </landmark>
    <landmark>
        <name>CNR-Green</name>
        <expected>58.220025 16.119286 -6.931881</expected>
        <detected>58.497813 9.654831 -7.491801</detected>
        <isPrecise>no</isPrecise>
        <fitResidual>6.319076</fitResidual>
    </landmark>
    <landmark>
        <name>10mm_0_01</name>
        <expected>-89.285621 28.732528 -32.248880</expected>
        <detected>-90.196775 29.579927 -32.380159</detected>
        <isPrecise>yes</isPrecise>
        <fitResidual>1.186873</fitResidual>
    </landmark>
    <landmark>
        <name>10mm_0_02</name>
        <expected>-87.561132 26.488978 27.684355</expected>
        <detected>-87.875622 26.872772 26.667286</detected>
        <isPrecise>yes</isPrecise>
        <fitResidual>0.919448</fitResidual>
    </landmark>

```

[more landmarks]

```

</landmarkList>
</phantom>

```

C Phantom Geometry

```

/*
 * Measurements were derived manually from the following document:
 * http://www.phantomlab.com/library/pdf/magphan\_adni\_manual.pdf
 * They can, therefore, be used without reference to ADNI publications.
 */
const cmtk::MagphanEMR051::SphereEntryType
cmtk::MagphanEMR051::SphereTable[cmtk::MagphanEMR051::NumberOfSpheres] =
{
    //
    // LICENSING EXCEPTION
    // Unlike the remainder of this file, the table of phantom sphere coordinates
    // is licensed under the CC BY 3.0 license
    // (https://creativecommons.org/licenses/by/3.0/us/)
    //
    // 1x 6.0cm SNR sphere
    { "SNR", 60, { 0.0, 0.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    // 2x 1.5cm spheres
    { "15mm@90mm", 15, { 89.0, -2.9, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "15mm@60mm", 15, { 0.0, -2.9, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    // 4x 3.0cm spheres -- estimated y coord's -- these are not marked in construction drawing
    { "CNR-Orange", 30, { 60.0, 20.0, 0 }, 0.590, 450, Self::SPHERE_COLOR_ORANGE },
    { "CNR-Red", 30, { -60.0, 20.0, 0 }, 0.430, 600, Self::SPHERE_COLOR_RED },
    { "CNR-Yellow", 30, { -60.0, -20.0, 0 }, 0.295, 750, Self::SPHERE_COLOR_YELLOW },
    { "CNR-Green", 30, { 60.0, -20.0, 0 }, 0.220, 900, Self::SPHERE_COLOR_GREEN },
    // 158x 1.0cm spheres
    // Plane 0
    // outer ring
    { "10mm_0_01", 10, { -86.4, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_02", 10, { -86.4, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_03", 10, { 86.4, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_04", 10, { 86.4, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_05", 10, { -64.7, 0.0, 64.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_06", 10, { -64.7, 0.0, -64.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_07", 10, { 64.7, 0.0, 64.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_08", 10, { 64.7, 0.0, -64.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_09", 10, { -30.0, 0.0, 86.4 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_10", 10, { -30.0, 0.0, -86.4 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_11", 10, { 30.0, 0.0, 86.4 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_12", 10, { 30.0, 0.0, -86.4 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_13", 10, { 0.0, 0.0, 91.5 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_14", 10, { 0.0, 0.0, -91.5 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    // middle ring
    { "10mm_0_15", 10, { -30.0, 0.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_16", 10, { -30.0, 0.0, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_17", 10, { 30.0, 0.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_18", 10, { 30.0, 0.0, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_19", 10, { -60.0, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_20", 10, { -60.0, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_21", 10, { 60.0, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    { "10mm_0_22", 10, { 60.0, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
    // inner ring

```

```

{ "10mm_0_23", 10, { -30.0, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_0_24", 10, { -30.0, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_0_25", 10, { 30.0, 0.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_0_26", 10, { 30.0, 0.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      single inferior mid-sagittal sphere
{ "10mm_0_27", 10, { 0.0, 0.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      single right sphere
{ "10mm_0_28", 10, { -91.5, 0.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 1
//      mid-sagittal
{ "10mm_1_01", 10, { 0.0, -30.0, 40.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_02", 10, { 0.0, -30.0, -40.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_03", 10, { 0.0, -30.0, 86.5 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_04", 10, { 0.0, -30.0, -86.5 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_05", 10, { 0.0, -30.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_06", 10, { 0.0, -30.0, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 30mm lateral
{ "10mm_1_07", 10, { 30.0, -30.0, 81.1 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_08", 10, { 30.0, -30.0, -81.1 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_09", 10, { -30.0, -30.0, 81.1 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_10", 10, { -30.0, -30.0, -81.1 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_11", 10, { 30.0, -30.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_12", 10, { 30.0, -30.0, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_13", 10, { -30.0, -30.0, 60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_14", 10, { -30.0, -30.0, -60.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_15", 10, { 30.0, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_16", 10, { 30.0, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_17", 10, { -30.0, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_18", 10, { -30.0, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 60mm lateral
{ "10mm_1_19", 10, { 60.0, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_20", 10, { 60.0, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_21", 10, { -60.0, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_22", 10, { -60.0, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 61mm lateral
{ "10mm_1_23", 10, { 61.0, -30.0, 61.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_24", 10, { 61.0, -30.0, -61.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_25", 10, { -61.0, -30.0, 61.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_26", 10, { -61.0, -30.0, -61.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 81.1mm lateral
{ "10mm_1_27", 10, { 81.1, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_28", 10, { 81.1, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_29", 10, { -81.1, -30.0, 30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_30", 10, { -81.1, -30.0, -30.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 86.5mm lateral
{ "10mm_1_31", 10, { 86.5, -30.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1_32", 10, { -86.5, -30.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 1b
//      +- 15mm lateral
{ "10mm_1b_01", 10, { 15.0, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_02", 10, { 15.0, 29.1, 65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_03", 10, { 15.0, 29.1, 85.2 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_04", 10, { 15.0, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_05", 10, { 15.0, 29.1, -65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },

```

```

{ "10mm_1b_06", 10, { 15.0, 29.1, -85.2 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_07", 10, { -15.0, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_08", 10, { -15.0, 29.1, 65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_09", 10, { -15.0, 29.1, 85.2 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_10", 10, { -15.0, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_11", 10, { -15.0, 29.1, -65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_12", 10, { -15.0, 29.1, -85.2 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 40mm lateral
{ "10mm_1b_13", 10, { 40.0, 29.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_14", 10, { 40.0, 29.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_15", 10, { -40.0, 29.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_16", 10, { -40.0, 29.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 45mm lateral
{ "10mm_1b_17", 10, { 45.0, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_18", 10, { 45.0, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_19", 10, { -45.0, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_20", 10, { -45.0, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_21", 10, { 45.0, 29.1, 65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_22", 10, { -45.0, 29.1, -65.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_23", 10, { 45.0, 29.1, -73.9 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_24", 10, { -45.0, 29.1, 73.9 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 64.5mm lateral
{ "10mm_1b_25", 10, { 64.5, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_26", 10, { -64.5, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 73.9mm lateral
{ "10mm_1b_27", 10, { 73.9, 29.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_28", 10, { -73.9, 29.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 85.2mm lateral
{ "10mm_1b_29", 10, { 85.2, 29.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_30", 10, { 85.2, 29.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_31", 10, { -85.2, 29.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_1b_32", 10, { -85.2, 29.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 2
//      +- 15mm lateral
{ "10mm_2_01", 10, { 15.0, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_02", 10, { 15.0, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_03", 10, { -15.0, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_04", 10, { -15.0, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_05", 10, { 15.0, -60.0, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_06", 10, { 15.0, -60.0, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_07", 10, { -15.0, -60.0, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_08", 10, { -15.0, -60.0, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_09", 10, { 15.0, -60.0, 67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_10", 10, { 15.0, -60.0, -67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_11", 10, { -15.0, -60.0, 67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_12", 10, { -15.0, -60.0, -67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 45mm lateral
{ "10mm_2_13", 10, { 45.0, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_14", 10, { 45.0, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_15", 10, { -45.0, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_16", 10, { -45.0, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 48.7mm lateral
{ "10mm_2_17", 10, { 48.7, -60.0, 48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_18", 10, { 48.7, -60.0, -48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },

```

```

{ "10mm_2_19", 10, { -48.7, -60.0, 48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_20", 10, { -48.7, -60.0, -48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 67.3mm lateral
{ "10mm_2_21", 10, { 67.3, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_22", 10, { 67.3, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_23", 10, { -67.3, -60.0, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2_24", 10, { -67.3, -60.0, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 2b (same as Plane 2 but at y=+59.1)
//      +- 15mm lateral
{ "10mm_2b_01", 10, { 15.0, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_02", 10, { 15.0, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_03", 10, { -15.0, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_04", 10, { -15.0, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_05", 10, { 15.0, 59.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_06", 10, { 15.0, 59.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_07", 10, { -15.0, 59.1, 45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_08", 10, { -15.0, 59.1, -45.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_09", 10, { 15.0, 59.1, 67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_10", 10, { 15.0, 59.1, -67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_11", 10, { -15.0, 59.1, 67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_12", 10, { -15.0, 59.1, -67.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 45mm lateral
{ "10mm_2b_13", 10, { 45.0, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_14", 10, { 45.0, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_15", 10, { -45.0, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_16", 10, { -45.0, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 48.7mm lateral
{ "10mm_2b_17", 10, { 48.7, 59.1, 48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_18", 10, { 48.7, 59.1, -48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_19", 10, { -48.7, 59.1, 48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_20", 10, { -48.7, 59.1, -48.7 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      +- 67.3mm lateral
{ "10mm_2b_21", 10, { 67.3, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_22", 10, { 67.3, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_23", 10, { -67.3, 59.1, 15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_2b_24", 10, { -67.3, 59.1, -15.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 3
{ "10mm_3_01", 10, { 28.3, -82.2, 28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_02", 10, { 28.3, -82.2, -28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_03", 10, { -28.3, -82.2, 28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_04", 10, { -28.3, -82.2, -28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_05", 10, { 0.0, -88.0, 25.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_06", 10, { 0.0, -88.0, -25.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_07", 10, { 25.0, -88.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_08", 10, { -25.0, -88.0, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3_09", 10, { 0.0, -89.5, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
//      Plane 3b
{ "10mm_3b_01", 10, { 28.3, 81.4, 28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_02", 10, { 28.3, 81.4, -28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_03", 10, { -28.3, 81.4, 28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_04", 10, { -28.3, 81.4, -28.3 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_05", 10, { 0.0, 87.2, 25.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_06", 10, { 0.0, 87.2, -25.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_07", 10, { 25.0, 87.2, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },

```

```
{ "10mm_3b_08", 10, { -25.0, 87.2, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE },
{ "10mm_3b_09", 10, { 0.0, 88.6, 0.0 }, 0.820, 282, Self::SPHERE_COLOR_NONE }
//
// END LICENSING EXCEPTION
//
```